

Total Service for Customer
The Leading Partner of Bioscience Technologies



We are searching for pureness,
overcoming limitness

EXL[®]7
Lab Water System



01 Time & Date

Year/Month/Day/Hour/Minute/Day of the week

02 Tap water quality indicator (Conductivity, Temperature)

By monitoring the supply status of the tap water, it secures the safety to protect the equipment in case of a shortage.

03 UV lamp light on/off

The real-time operation status of the UV lamp maintains optimum condition.

04 The name & The life of filter

It shows filter use, history and details.

05 Start & Stop button Pure/Ultrapur Water dispensing

06 Pure Water status indicator

Displays Conductivity, Temperature, Rejection rate (%)
System controlled operation has voice guidance and pop-up message whenever it is out of range set point under monitoring real time.

Ultrapure Water status indicator

Displays Conductivity, Temperature, TOC (Option)
System controlled operation has voice guidance and pop-up message whenever it is out of range set point under monitoring real time.

07 Tank level indicator

It is easy to check tank level displayed in real time.

08 Pure /Ultrapure Water producing volume

09 Pure /Ultrapure Water producing volume control button

User can produce off water as wanted. (0.1L ~99.9L)

10 Menu (System information indicator)

Displays Model, S/N, Software version, System alarm history, Setting.

11 Print out button

System can print out data.

(ex. Used & Available time of filter, Filter usage, Filter exchange history, Alarm history, Water quality)

12 Pre filter flushing button

13 Pure /Ultrapure Water producing button

14 Pure /Ultrapure Water producing indicator

Auto RO Flushing

It can flush RO filter automatically to control optimal filter quality before and after producing Pure Water.

History storage

Consumable exchange, Alarm, Set point

Data storage

Filter Information (Used & Available time, Total output), Alarm history, Consumable exchange history are saved in USB memory.

TOC (Option)

Displays TOC value of Ultrapure Water and a message pops up on screen when it is out of set point.

Protection for Ultrapure Water filter and pump

System allows Ultrapure Water filter to stop automatically and alarm when water level in Pure Water storage tank falls below.

Software update

It is easy to be up to date software by USB memory anybody.

Remote control

Direct system control via internet network & technician's direct access.

POU (Point Of Use Gun) option

Built in quantitative sampling function to produce off Pure and Ultrapure Water as much as user wants. (0.1L ~ 99.9L units)



EXL[®]7

We are searching for pureness, overcoming limitness

Vivagen Co., Ltd. developed EXL[®]7 Lab Water System which consists of experts who have more than 10 years of experience in R&D. We are sure our expertise and technical know-how help customers solve their complex problems and achieve their goals.

Pure/Ultrapure Water System

EXL[®]7 can be used for both Pure Water and Ultrapure Water simultaneously.

In particular, Pure Water can be produced while dispensing Ultrapure Water, which allows for convenient experiment.

Fast Flow Rate

Pure Water is up to 16L(32L)/hour and Ultrapure Water is up to 120 L/hour.

Touch screen (7") & Voice guidance

For ease of use, the system is installed touch screen, the Setting of voice and Pop-up message guidance.

Built-In 20L Pressurized Tank

While open tanks from other companies have the risk of contamination by air, EXL[®]7 uses pressurized tank which can perfectly block external contamination. Beyond basic inner tank, you can have option available outer tank 20L, 30L, 70L without additional sensor, the height of tank.

High Capacity DI Cartridge

It uses the filter with excellent processing capacity, which can assure best water quality.

UV Lamp & Leakage detector installed by default

EXL[®]7 has a UV lamp installed by default, which can produce Ultrapure Water with extremely low level of organic substances.

A leakage detector is installed by default, which issues alarming sound and message if water leakage occurs to ensure user safety.


Touch screen monitor

- System information ● The life of filter ● UV lamp indicator ● Tap water quality indicator
- The status of Pure /Ultrapure Water dispensing ● Pump monitor ● Pure /Ultrapure Water producing volume
- Total amount ● Print out ● Simple installation & initial setting guide ● Auto Flushing ● Record history
- USB port ● TOC(Optional) ● Software upgrade ● Tank level indicator ● Remote control

Foot switch (Option)

The system provides foot switch against busy two hands.



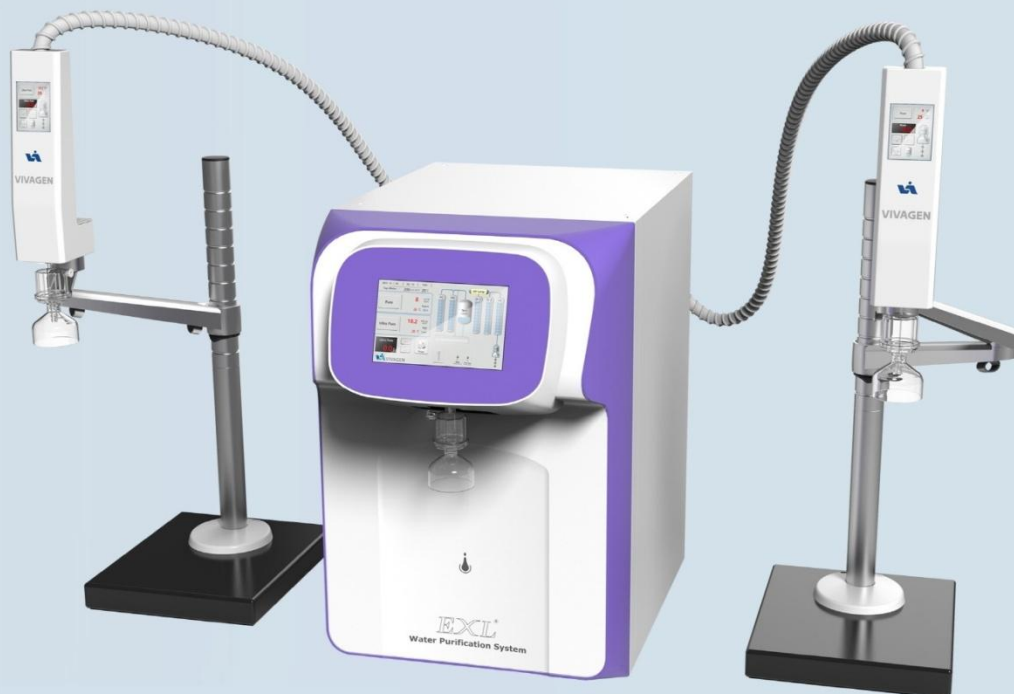

EXL®7 Bio

ASTM Type III&I Grade (Bio)

• Ultra-trace and trace-organic and organic analysis • HPLC • GC-MS • ICP-MS • TOC analysis
• Ion chromatography • Solid phase extraction • Electrochemistry

• Molecular biology • PCR • DNA sequencing • Cell culture • Media preparation • Microbiology
• In-vitro fertilization • Electrophoresis • Monoclonal antibody production

EXL®7 Lab Water System Pure & Ultrapure Water Quality

Series	Resistivity (at 25°C)	• TOC	Ultrapure						RO Module Rejection Rates (Pure)							
			•• Bacteria	••• Endotoxin	••• RNases (RNase-free)	••• DNases (DNase-free)	Particles (@0.2µs)	Flow Rate	Resistivity (at 25°C)	Monovalent Ions	Polyvalent Ions	Particles	Microorga- nisms	Dissolved Organics (>300MW)	Product Flow Rate (±20%at25°C)	Maximum Recovery Rate*
	MΩxcm	ppb	CFU/mL	EU/mL	pg/mL	pg/mL	/ml	(l/min)	MΩ.cm	%	%	%	%	%	l/h	%
EXL®7 Bio	18.2	<5	<0.01	<0.001	<1	<5	<1	Up to 2	Up to 98	Up to 98	>99	>99	>99	>99	Up to 16/32	Approx. 50

Depends on feed water quality and temperature.

* TOC specs - Test Conditions: EXL®System equipped with Pre Carbon filter and DI cartridge and with feed (tap) water quality within specifications.
Product water quality may vary due to local feed water conditions.

** Results with 0.2 Final Filter.

*** Results with Ultra Filter polisher in place.

Feed Water Quality

*Exclusively potable water according to the drinking water regulations of the USA, the European Union, Korea or Japan.

Temperature	Conductivity (at 25°C)	TDS	TOC	Max. Permanent Hardness (max. CaCO ₃)	Free chlorine	Iron (Total Fe Content)	Fouling Index	Turbidity	pH Value	Input Pressure
°C	µS/cm	ppm	ppb	ppm	ppm	ppm	SDI	NTU		
2 – 30	< 500	< 250	< 2000	360	3	< 0.1	< 5	< 1	4 – 10	0.5 – 6.9 bar, Recommended > 2 bar

Specifications

	EXL®7 Bio
Power	220~240VAC/ 50/60Hz (100W)
System Weight	37kg
System Dimension	375 (W) X 600 (D) X 600 (H) mm
Packing Dimension	500 (W) X 700 (D) X 800 (H) mm
Operation Temperature & Humidity Range	2~35 °C / Under 80% Relative Humidity
Storage Temperature & Humidity Range	5~35 °C / Under 80% Relative Humidity

Founded in 2003, Vivagen Co., Ltd. is the company which consists of experts who have more than 10 years of experience in sales, marketing and R&D.

Our products consist of self-developed Pure & Ultrapure integrated water system (EXL® Series), Gas control module for cell culture (EXL® Navigator), Transfection Reagent (VivaMagic), Taq Polymerase and Marker. The patent of Lab Water System is registered with Korean Intellectual Property Office in Oct. 2010. In 2016, we released new product, EXL®7.

Vivagen Co., Ltd. provides differentiated services to customers in various experiments such as filtration, purification, concentration, incubation and analysis.

We are sure our expertise and technical know-how help customers solve their complex problems and achieve their goals.



Ordering Information

Order No.	Description	Type
System		
EXL®7 Bio 16	EXL®7 Bio 16 Pure & Ultrapure Water System(w/UV lamp) & 20L Tank	III & I Grade
EXL®7 Bio 32	EXL®7 Bio 32 Pure & Ultrapure Water System(w/UV lamp) & 20L Tank	III & I Grade
Option		
LM-VAL T	IQ/OQ Validation EXL®7 Pure & Ultra Pure	
LM-EXL7-ITOC	Inline TOC for EXL 7	
LM-EXL-Printer	Printer for EXL5&7	
LM-EXL-POU	Piont of use gun Pure or Ultra Pure for EXL5&7	
LM-EXL-FW	Foot Switch EXL5&7	
LM-20Tank	20L Pressure Tank	
LM-30Tank	30L Pressure Tank	
LM-70Tank	70L Pressure Tank	
Consumables		
LM-EXL7-PC	Pre Carbon Filter for EXL7 (4ea/pk)	
LM-EXL7-RO	RO Membrane for EXL7	
LM-EXL7-DI	DI Cartridge#1 & #2 for EXL7 (each 2ea/pk)	
LM-EXL-UV	UV Lamp for EXL3&5&7	
LM-EXL7-UF	UF Membrane for EXL7	
LM-FINAL	0.2u Final filter for EXL3&5&7	
LM-EXL7-1M	1um Pre filter for EXL7 (4ea/pk)	
LM-EXL7-5M	5um Pre filter for EXL7 (4ea/pk)	



Suite 618, Seongnam Woorim Lions Valley I, 27, Dunchon-daero 457beon-gil, Jungwon-gu, Seongnam-si, Gyeonggi-do, Republic of Korea, 13219
Tel) +82-31-737-2080 Fax) +82-31-737-2083 / Visit us at www.VIVAGEN.co.kr
E-mail : vivagen@vivagen.co.kr / Shopping Mall: www.vgmall.co.kr

2018 VIVAGEN All right reserved. This product or document is protected by copyright and distributed under licenses restricting its use, copying, distribution, and decompilation. No part of this product or document may be reproduced in any form by any means without prior written authorization of VIVAGEN and its licensors, if any. Products shown on this catalog are subject to change without any prior notice.